


## Method

Method ID	00002r1
Internal ID	M55
Valid	Yes
User defined label 1	
Can run on blocked instr.	No

Name	Silbiger TA (EP)r1
Version	8
Release state	Released

Comment

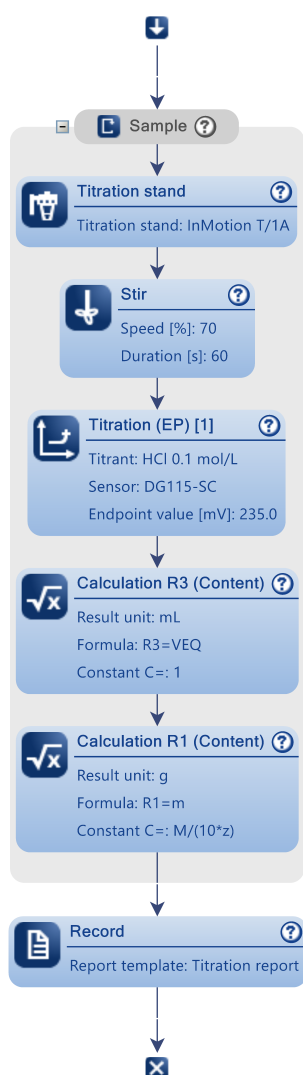

## Export templates used

Hawaii\_Samples\_1010201 CSV Simple  
9

## Sample (S1)

### Sample

Sample type	Sample
Number of IDs	1
ID 1	--
Entry type	Weight
Lower limit [g]	59
Upper limit [g]	61
Density [g/mL]	1.0
Correction factor	0
Temperature [°C]	25.0
Entry	Arbitrary
InMotion reader	None
Number of sample factors	0

## Titration stand (Titrationstand1)

### Titration Stand

Type	InMotion T/Tower A
Titration stand	InMotion T/1A
Head position	Sample
Lid handling	No

## Stir (Stir1)

### Stir

Speed [%]	70
Duration [s]	60
Condition	No

## Titration (EP) [1] (TitrationEP1)

### Titrant

Titrant	HCl
---------	-----

Concentration [mol/L]	0.1
-----------------------	-----

## Sensor

Type	pH
Sensor	DG115-SC
Unit	mV

## Temperature Acquisition

Temperature Acquisition	Yes
Temperature sensor	DT1000
Temperature unit	°C

## Stir

Speed [%]	30
-----------	----

## Predispense

Mode	Potential
Potential [mV]	185
Wait time [s]	240

## Control

Endpoint type	Absolute
Tendency	None
Endpoint value [mV]	235.0
Control band [mV]	50
Dosing rate (max) [mL/min]	0.5
Dosing rate (min) [μL/min]	10

## Termination

At EP	Yes
Termination delay [s]	0
At Vmax [mL]	4.0
Max. time infinite	Yes
Max. time [s]	∞

## Accompanying Stating

Accompanying stating	No
----------------------	----

## Condition

Condition	No
-----------	----

## Content (R3)

### Calculation

Result	Content
Result unit	mL
Formula	R3=VEQ
Constant C=	1
M	M[Hydrochloric acid]
z	z[Hydrochloric acid]

Decimal places	3
Result limits	No
Extra statistical functions	No
Send to buffer	No
Write to Smart Tag	None
Condition	No

## Content (R1)

### Calculation

Result	Content
Result unit	g
Formula	$R1=m$
Constant C=	$M/(10*z)$
M	M[Hydrochloric acid]
z	z[Hydrochloric acid]
Decimal places	3
Result limits	No
Extra statistical functions	No
Send to buffer	No
Write to Smart Tag	None
Condition	No

## Record (Record1)

### Report

Report template	Titration report
Print	No
Condition	No

## Signing History

Time	Full name	User login name	Role name	Operation	Comment