


## Method

Method ID	PH CAL
Internal ID	M22
Valid	Yes
User defined label 1	
Can run on blocked instr.	No

Name	Calibration
Version	5
Release state	Released

Comment

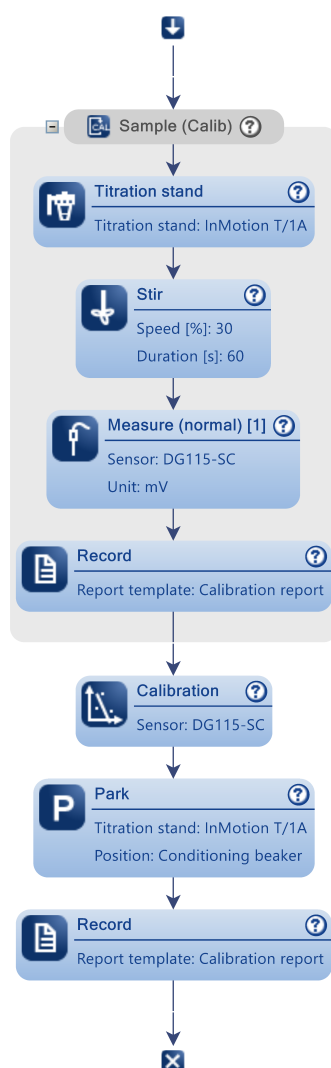

## Export templates used

## Sample (Calib) (C1)

### Sensor

Sensor type	pH
Sensor	DG115-SC
Unit	pH

### Calibration

Action	Calibration
Automatic buffer recognition	No
pH buffer list	METTLER TOLEDO USA (Ref. 25°C)
Calibration	Linear
Number of buffers	3
Buffer 1	4.01
Buffer 2	7.00
Buffer 3	10.01

### Temperature

Temperature [°C]	25.0
------------------	------

## Titration stand (Titrationstand1)

### Titration Stand

Type	InMotion T/Tower A
Titration stand	InMotion T/1A
Head position	Sample
Lid handling	No

## Stir (Stir1)

### Stir

Speed [%]	30
Duration [s]	60
Condition	No

## Measure (normal) [1] (Measurenormal1)

## Sensor

Type	pH
Sensor	DG115-SC
Unit	mV

## Temperature Acquisition

Temperature Acquisition	No
-------------------------	----

## Stir

Speed [%]	30
-----------	----

## Acquisition of Measured Values

Acquisition	Equilibrium controlled
dE [mV]	0.5
dt [s]	1
t(min) [s]	10
t(max) [s]	30
Mean value	No

## Condition

Condition	No
-----------	----

## Record (Record1)

### Report

Report template	Calibration report
Print	No
Condition	No

## Calibration (Calibration1)

### Calibration

Sensor	DG115-SC
Limits	No
Condition	No

## Park (Park1)

### Park

Titration stand	InMotion T/1A
Position	Conditioning beaker
Condition	No

## Record (Record2)

### Report

Report template	Calibration report
Print	No
Condition	No

## Signing History

Time	Full name	User login name	Role name	Operation	Comment